

pluspoint FAST

Functional Knee orthosis



**ORTHOSERVICE
RO+TEN**

Take care feel better

The all-time best support for the knee

Designed to support knee instability or torn ligaments, it represents a first choice brace for **any instability of the knee**, whether it is due to traumatic or degenerative reasons. Recommended for both **preop and postop phases**, as a conservative treatment or as **preventive protection**.

Front fit

The Pluspoint Fast knee brace is put on from the front without having to insert the foot. This feature makes it extremely **practical and easy to use**.

THE STRAPS ARE INCREDIBLY EASY TO FASTEN

The straps don't need to be removed and reinserted every time, because they fasten with **small levers**, which maintain the regulation made on the first application.




EFFICIENCY OF THE 4 THRUST POINTS

The 4 thrust points system, guarantees efficient support, providing optimal stabilization of the knee.

pluspoint FAST

Shapable anatomical frame with front fit

Made of a **special aluminum alloy**, it is shapable, especially on the shin and thigh. The **front fit** makes it easy for the patient: fast and intuitive to wear.

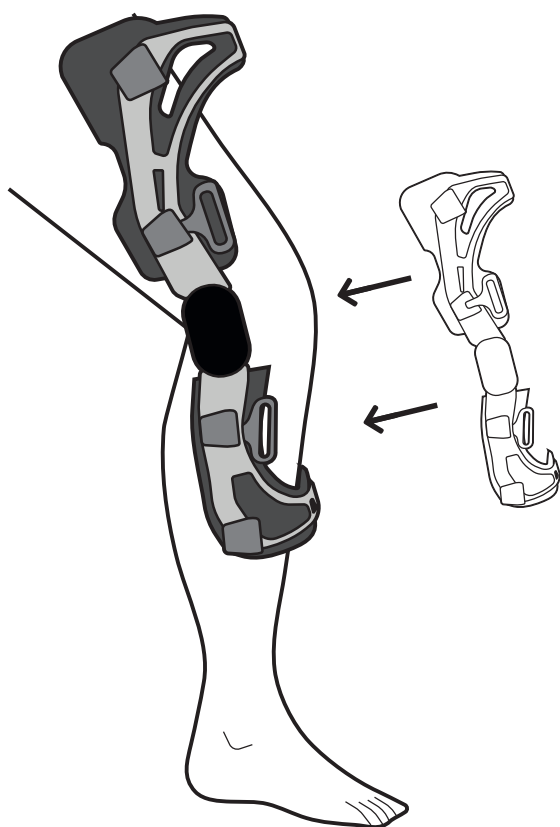



CLOSURE WITH LEVER BUCKLES

The straps are not inserted in the buckles but are hooked to the frame by means of **special lever buckles**, also equipped with **safety buttons**.



NUMBERED AND PADDED STRAPS

The straps are **numbered** to facilitate the correct closing sequence. Internally they have **MTP padding** and are **seamless**, for more comfort on the skin.

MEDIAL CONDYLAR PAD

Anatomically designed and pre-assembled, to provide greater stability.

SET OF CONDYLAR PADS

A practical set of pads, of **different thicknesses**, is included, suitable for any need.

TIBIAL PADDING

The tibia is positioned between two additional lateral shims, to enhance **the comfort and the stability** of the brace.

Never before
has a knee brace
been so easy to wear

CLICKCLACK JOINT

It adjusts the flexion-extension by means of wedges, **without the need of any tool** (screwdriver or allen wrench). Polycentric and with a very low profile, you open it, adjust it, and close it.

Extension: 0°, 10°, 20°, 30°, 40°,

Flexion 45°, 60°, 75°, 90°

Blocks at 0°, 10°, 20°, 30°.



THERMOFORMED PADS

Internally (where in contact with the skin) they are covered with a hypoallergenic and non-slip fabric, which has a modern texture design **externally**. **Fastens** with low-profile snaps that are easily removed to clean the brace.



pluspoint FAST

Functional Knee orthosis

REF. 22040

Indications

- ACL, MCL, LCL, PCL instability and ruptures, combined instability
- Complex instabilities of the knee joint, traumatic or degenerative
- Pre/Post-surgery both rehabilitative or conservative
- Ligament protection during sports

SINGLE PATIENT
MULTIPLE USEPATENT
EP3492050PATENT
US10548758EASY
ADJUSTINGEASY
FASTENINGLOW
PROFILEANATOMICAL
FITCUSTOMIZABLE
FIT

Item

Size

Thigh circ. taken 15 cm above
mid-patella

Length cm

Colour

REF.22040

S

39/47

M

47/53

L

53/60

XL

60/67

XXL

67/75

40

40

40

40

40

dark grey

Specify right or left



orthoservice.com



Headquarter: ORTHOSERVICE AG

Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland

Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89

info@orthoservice.com · www.orthoservice.com

Niederlassung Deutschland: Orthoservice Deutschland GmbH

Flugstraße 8 · D-76532 Baden-Baden · Deutschland

Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13

info@orthoservice.de · www.orthoservice.de

Siedziba w Polsce: Orthoservice Polska Sp. z o. o.

ul. Warszawska 416a · 42-209 Częstochowa · Polska

Tel. +48 (0) 34 340 13 10

info@orthoservice.pl · www.orthoservice.pl

Sede italiana: RO+TEN S.r.l.

Sede legale: Via Marco De Marchi, 7 · I-20121 Milano (MI) · Italia

Sede operativa e amministrativa:

Via Comasina, 111 · I-20843 Verano Brianza (MB) · Italia

Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34

info@roplusten.com · www.roplusten.com

Società soggetta a Direzione e Coordinamento (art. 2497bis CC):

Orthoservice AG (CH) · 6830 Chiasso (TI) · Switzerland


**ORTHOSERVICE
RO+TEN**

Take care feel better