

9000P00196INT01.0324

elbo PRO

M 24017R
M 24017L

Telescopic elbow brace with graduated movement

Code	M 24017R · M 24017L
Size	universal
Length of the hinges cm	36-46
Colour	dark grey

Right or left




Indications

- Compound fractures of elbow, distal humerus, proximal radius and ulna
- Post-operative immobilisation following tendon and ligament repairs
- Useful in the rehabilitation phase following surgical and conservative treatment of elbow fractures
- Elbow subluxations and dislocations
- Chronic elbow injuries
- Elbow hyperextension
- Post-traumatic treatment



elbo PRO

Telescopic elbow brace with graduated movement

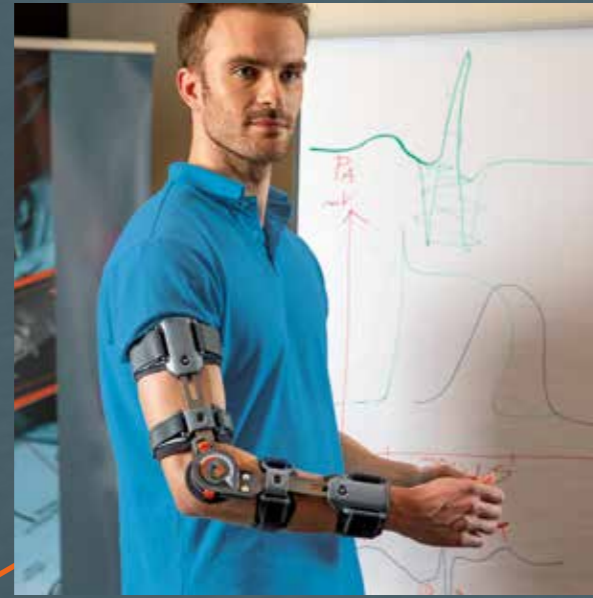

...satisfies all sizes in case of:

- Compound fractures of the elbow, distal humerus, proximal radius and ulna
- post-operative immobilisation following tendon and ligament repairs
- following surgical and conservative treatment of elbow fractures
- elbow subluxations and dislocations
- chronic elbow injuries
- elbow hyperextension
- post-traumatic treatment



orthoservice.com



Headquarter: ORTHOSERVICE AG
 Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
 Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
 info@orthoservice.com · www.orthoservice.com

ORTHO SERVICE DEUTSCHLAND GmbH
 Flugstraße 8 · D-76532 Baden-Baden · Deutschland
 Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13
 info@orthoservice.de · www.orthoservice.de

ORTHO SERVICE POLSKA Sp. z o.o.
 ul. Warszawska 416a · 42-209 Częstochowa · Polska
 Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785
 info@orthoservice.pl · www.orthoservice.pl

RO+TEN S.r.l.
 Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia
 Sede operativa e amministrativa:
 Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia
 Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34
 info@roplusten.com · www.roplusten.com



ORTHO SERVICE RO+TEN
 Take care feel better



ORTHO SERVICE RO+TEN
 Take care feel better

Efficient and versatile...

TELESCOPIC

Easily adapts to every arm eliminating the need for different sizes

THE JOINT IS EASY AND SAFE

It is regulated without the need for tools but cannot be unintentionally released




TELESCOPIC LENGTH ADJUSTMENT OF THE ARTICULATED HINGE

thanks to easy-to-use buttons on the valves



FIRM ANCHORAGE

(the brace doesn't rotate) thanks to the polyamide valves with mouldable aluminium arm grip

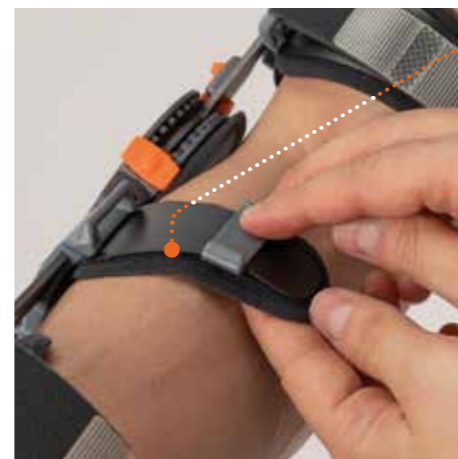

MAXIMUM COMFORT

thanks to the heat-sealed padding in soft lycra that comes into contact with the skin



EASY FLEXION-EXTENSION ADJUSTMENT

thanks to the 'Pull&Set' joint with spring-loaded sliders (pull outwards, position to desired degree, release)

Elbo PRO, the new era for elbow braces

Included accessories

LENGTH-ADJUSTABLE SHOULDER STRAP

with padded shoulder protection to be worn on the ipsilateral shoulder for optimal comfort



NEUTRAL HAND SUPPORT

for forearm rotation control and a feeling of stability

THE SET FLEXION-EXTENSION

cannot be altered by arm movements or impacts against surfaces