

genuFEEL OA

**Gonarthrosis
knee brace**



**ORTHOSERVICE
RO+TEN**

Take care feel better

Time, a key factor

Relieving pain improves patients' quality of life:

- inoperable
- awaiting surgery
- at an age not yet suitable for surgical treatment


this is the real goal.

GenuFeel OA alleviates the side effects of degenerative gonarthrosis.



Fastening, hooking the straps, **lowering** the levers: a sequence of three simple gestures, made possible by innovative design, sophisticated materials and increasingly functional details.

The only patented system that takes less than 1 minute...

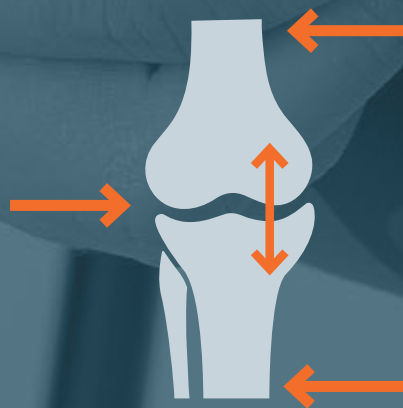


Knee unloading mechanism



...and promotes the rehabilitation of the entire joint

The real benefit of a brace working on **3 thrust points** is its deep action: the **thrust on the knee**, combined with the **action of the 4 straps** and the **contralateral lever system**, makes the joint assume a 'position' closer to its physiological one:

- **Immediate relief from local pain and relief of joint tendons**
- **A more natural posture**
- **Balanced walking**
- **to the benefit of the entire musculoskeletal system**

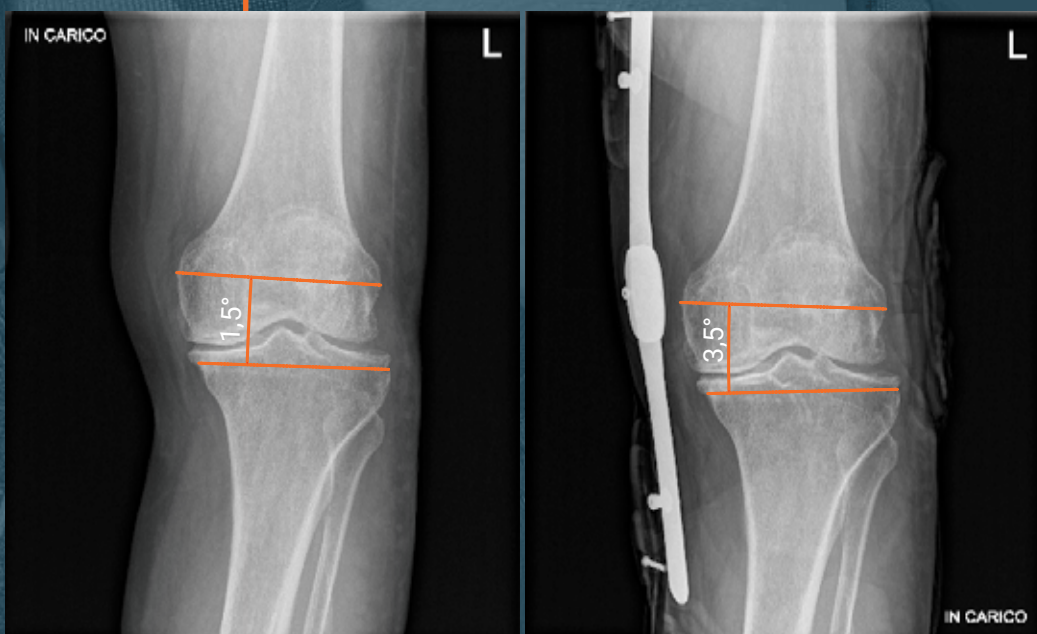

The comparative Rx data, with and without the brace worn, shows a varus correction of a maximum of 2 degrees, although numerically the data is to be considered small, it is nevertheless significant for all cases of inoperability or of a patient who is not ready for surgery (refusal of surgery).

Considering the clinical case of a major gonarthrosis and therefore, at least radiographically, already of surgical relevance, the correction of 2 degrees is to be considered valid in order to obtain the patient's clinical satisfaction by relieving the most important symptom: pain.

TOTAL COMFORT

on the skin, thanks to the **stretchy, breathable fabric**, made with no seams and no gluing, and to the padding incorporated in the pressure zones




EFFECTIVE DISTRIBUTION OF THRUST PRESSURE

thanks to the new ring which also features an internal protective padding, in **gel** and removable



FLEXION-EXTENSION ADJUSTMENT

with simplified joint

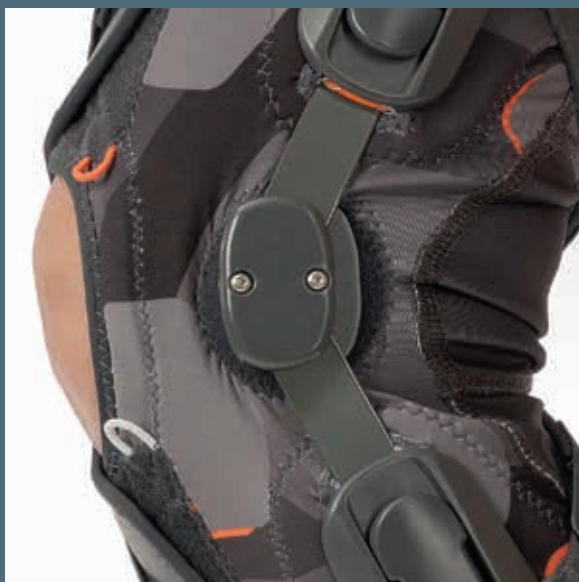



EASY TENSIONING OF THE STRAPS

thanks to **low-profile levers**
that only need a light pressure
from the palm of the hand

IMMEDIATE DETERMINATION OF STRAP TENSION,

for possible adjustments:
graduated scale



EASY AND QUICK FITTING

front opening and finger pockets that facilitate the use of the Velcro® fasteners



INTUITIVE PRE-POSITIONING AND POSITIONING

of the straps: coloured loops



The gonarthrosis
knee brace that works
when you need it



**SECURE FIT TO
THE LIMB
AT ALL TIMES,**

thanks to the anatomical
shape, the technopolymer
stabilising **valves**, the
steel **rod** and the flexible
polymer **splints**

genuFEEL OA

Gonarthrosis knee brace

M24014 R

M24014 L

Indications

- Moderate unicompartmental gonarthrosis
- Post-operative cartilage or meniscus surgery

Code	M24014 R • M24014 L						
Size	XS	S	M	L	XL	XXL	3XL
Circumf. 15cm above knee center cm	33/39	39/45	45/52	52/58	58/66	67/74	74/81
Height cm	40	40	40	40	40	40	40
Colour	dark gray						


Right or left

RIGHT = for medial gonarthrosis of the right knee (varus right knee) or for lateral gonarthrosis of the left knee (valgus left knee)

LEFT = for medial gonarthrosis of the left knee (varus left knee) or for lateral gonarthrosis of the right knee (valgus right knee)

ANATOMICAL
FIT

BREATHABLE

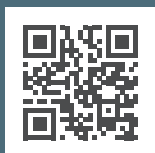LOW
PROFILEEASY
FASTENINGEASY
ADJUSTINGEASY
PULLINGCUSTOMIZABLE
FITSINGLE PATIENT
MULTIPLE USESOLO RIVENDITORI
AUTORIZZATI

www.orthoservice.com

SISTEMA DI GESTIONE
CERTIFICATO

Headquarter: ORTHOSERVICE AG
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
info@orthoservice.com · www.orthoservice.com



ORTHOSERVICE DEUTSCHLAND GmbH
Flugstraße 8 · D-76532 Baden-Baden · Deutschland
Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13
info@orthoservice.de · www.orthoservice.de

ORTHOSERVICE POLSKA Sp. z o.o.
ul. Warszawska 416a · 42-209 Częstochowa · Polska
Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785
info@orthoservice.pl · www.orthoservice.pl

OSUK Ltd

Units 47&50 North East BIC
Wearfield
Sunderland Enterprise Park East
Sunderland SR5 2TA
Tel +44 (0) 191 51 66 220
info@orthoservice.uk · www.orthoservice.uk

RO+TEN S.r.l.

Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia
Sede operativa e amministrativa:
Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia
Tel. 0039 0362 18 51 200
info@roplusten.com · www.roplusten.com



**ORTHOSERVICE
RO+TEN**

Take care feel better