

genupolar PRO

Stabilization,
compression and
cooling in a single
brace




ORTHOSERVICE
RO+TEN

Take care feel better

How it works

Three **synergistic actions** for optimal recovery

- 1 Stabilize:** the polycentric hinges help maintain the natural physiological alignment of the joint.
- 2 Cools:** the gel packs help reduce inflammation and pain.
- 3 Compresses:** the air chamber allows you to adjust the pressure, providing a draining and antiedema effect

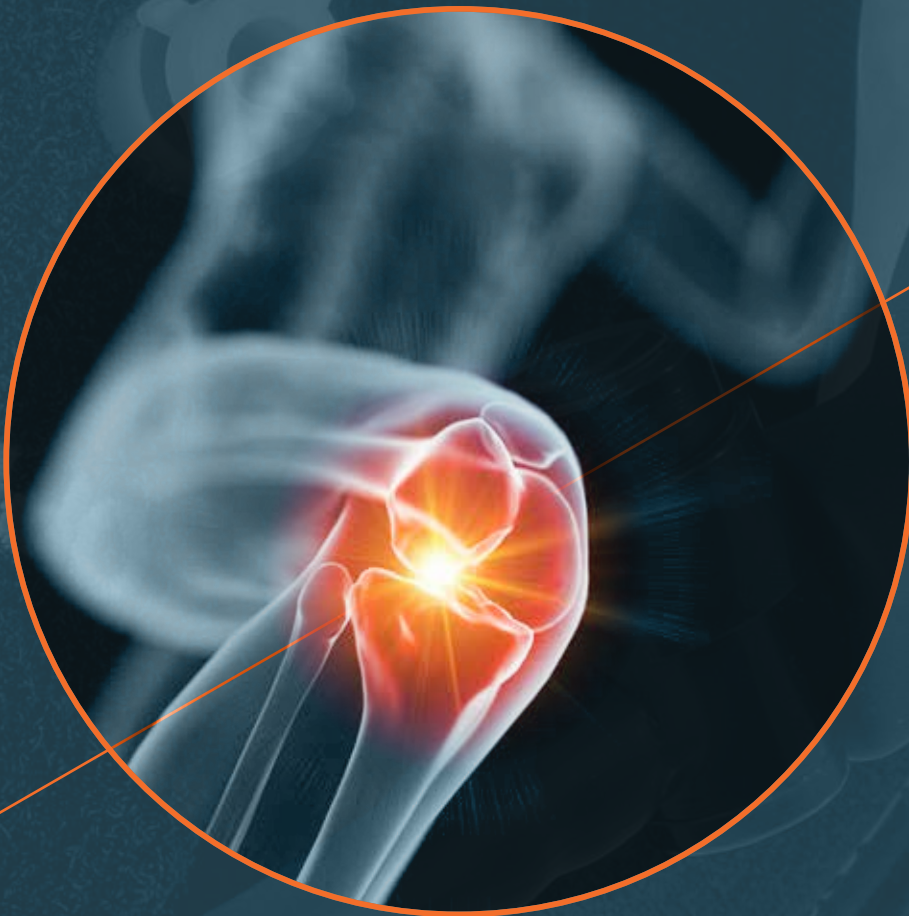

Why choose this brace

- Technology, comfort, and control for **optimal joint recovery**
- **Advanced solution** for knee treatment and rehabilitation
- **Modular** and **universal** system
- **Intuitive flexion–extension** adjustment
- **Easy** to apply and adjust

genupolar PRO

Mechanical stability and physical therapy to reduce swelling, pain, and edema



**SMART DESIGN
AND NEXTGEN
MATERIALS**

ONE SIZE FITS ALL

thanks to two dedicated
strap extensions
(included)

Cooling and compression accessories included to promote healing by effectively reducing swelling, edema, and pain



TWO REUSABLE GEL PACKS

Anatomically shaped with a patellar opening for cooling



INTEGRATED INFLATABLE AIR CHAMBER

with dedicated pump for personalized compression



EXTENSION STRAPS

for a circumference up to 75 cm




SOFT COMPRESSION SOCKS

The knee brace that speeds up the healing process

POLYCENTRIC HINGE

easily adjustable from the outside
using pins:

- Flexion 0°, 10°, 20°, 40°, 60°, 90°, 150°
- Extension 0°, 10°, 20°, 40°, 60°, 90°

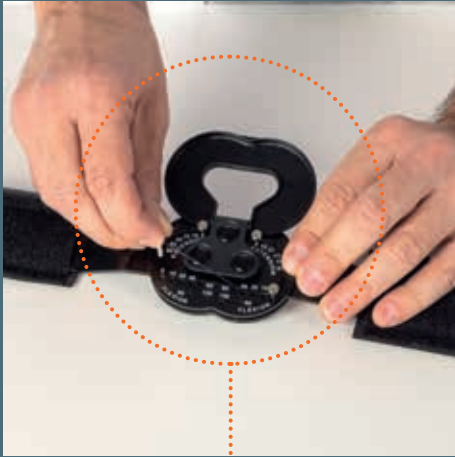

SOFT, DOUBLE LINED FOAM STRUCTURE

with integrated closure straps
to wrap around the knee
featuring a practical **Velcro®**
fastening

- **Outer layer**
velcrocompatible velour
- **Inner layer**
elastic fabric

genupolar PRO

REF. 26100

Knee brace with polycentric hinged rods for F-E control
and accessories for compression and cold therapy

Code	REF. 26100
Size	Universal
Length cm	34
Color	black

Fits right and left

Indications

Indicated in all cases requiring physiological stabilization of the knee joint, with the added possibility of cooling and compression (to help reduce inflammation, swelling, edema, and pain), such as:

- Mild to moderate collateral ligament instability
- Meniscal injuries
- Postsurgical recovery
- Osteoarthritis
- Patellofemoral pain

Contraindications

- DO NOT USE this product on patients with the following: vascular disorders, diabetes, sickle cell anemia, local circulatory disorders, tissue inflammation, hypersensitivity to cold, or previous cold injuries.



SINGLE PATIENT
MULTIPLE USE



Headquarter: ORTHOSERVICE AG
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
info@orthoservice.com · www.orthoservice.com



ORTHOERVICE DEUTSCHLAND GmbH
Flugstraße 8 · D-76532 Baden-Baden · Deutschland
Tel. +49 (0) 7221 40 533 0 · Fax +49 (0) 7221 40 533 20
info@orthoservice.de · www.orthoservice.de

ORTHOERVICE POLSKA Sp. z o.o.
ul. Warszawska 416a · 42-209 Częstochowa · Polska
Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785
info@orthoservice.pl · www.orthoservice.pl

OSUK Ltd
Units 47&50 North East BIC
Wearfield
Sunderland Enterprise Park East
Sunderland SR5 2TA
Tel +44 (0) 191 51 66 220
info@orthoservice.uk · www.orthoservice.uk

RO+TEN S.r.l.
Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia
Sede operativa e amministrativa:
Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia
Tel. 0039 0362 18 51 200
info@roplusten.com · www.roplusten.com



orthoservice.com



**ORTHOERVICE
RO+TEN**

Take care feel better