

# spinalplus 3

The spinal orthosis for osteoporosis enters in a new era



**ORTHOSERVICE  
RO+TEN**

Take care feel better



## Unique, patented tensioning system

### Pulling the terminals:

- all the posterior straps are tensioned
- the shoulders are adjusted simultaneously

### Practical benefits:

- faster adjustment
- uniform tension
- fewer set-up errors



### Maximum effectiveness with minimum effort:

- **multiplies** the force, **doubles** the traction on the shoulders
- **improves** posture control
- **increases** patient comfort


## Why choose Spinalplus 3

- 1 Telescopic and shapeable frame
- 2 Fast, centralised adjustment
- 3 Patented strap system that reduces manual effort
- 4 Guaranteed **stability** and practicality

## How does it work?

Spinalplus 3 is essential in the conservative/rehabilitative approach for **patients with osteoporosis or hyperkyphosis**.



The **rigid frame** and the **two counter-pressure points** at the shoulders and lower abdomen create a lever that supports the **thoracic spine**, limits flexion-extension, and helps straighten it.

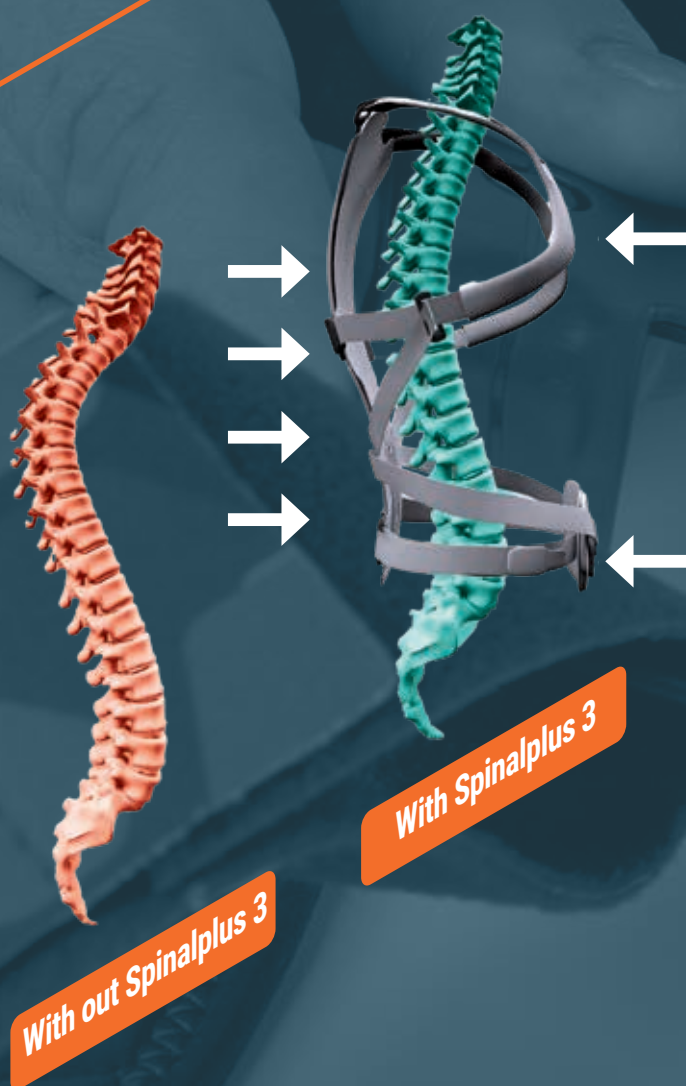

# Personalisation, lightness and adaptability

## TELESCOPIC FRAME

It can be easily lengthened or shortened through a **sliding system** and locking pin; it can be **shaped** according to the patient's anatomy

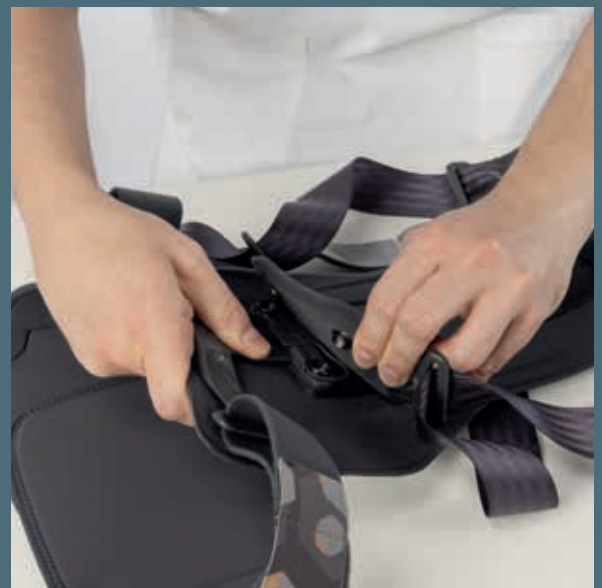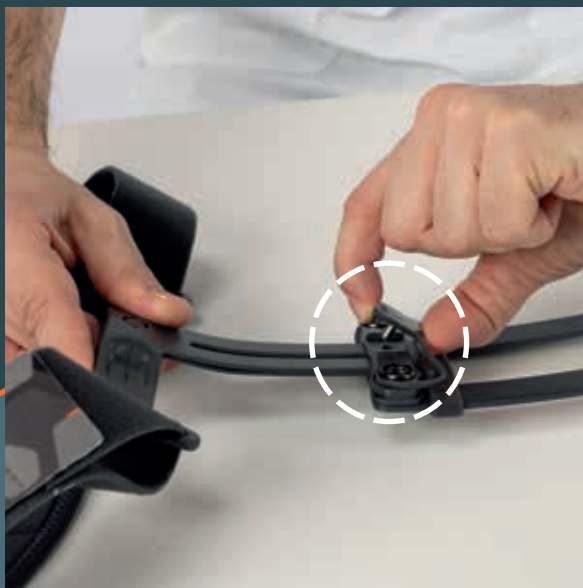

It can be removed, adjusted and reapplied in easy way thanks to **two practical locking buttons**



## SHOULDER STRAPS

with Velcro®-positionable padding

## SIMPLIFIED ADJUSTMENT

even when the frame is removed, the **straps** remain in position

## DORSAL PADDING

in bi-laminated PE

## EASY ADAPTATION

**Velcro®-adjustable** according to shoulder circumference

## ONLY 2 SIZES

to meet the needs of all patients

# Effective, anatomical and easy to wear



## **EASY TO CLOSE**

thanks to **Velcro®** system  
and integrated **handles**



## **A SINGLE ACTION**

complete control of all **terminals**



## **ABDOMINAL STABILITY**

thanks to semi-corset

Spinalplus 3:  
the evolution  
of the brace for  
osteoporosis



**COMFORTABLE,  
MODERN AND YET  
DISCREET**

it disappears under clothing

# spinalplus 3

## Spinal orthosis for osteoporosis

|                            |            |          |
|----------------------------|------------|----------|
| <b>Code</b>                | REF. 26049 |          |
| <b>Size</b>                | <b>1</b>   | <b>2</b> |
| <b>Distance T1-L5 (cm)</b> | 35/47      | 47/61    |
| <b>Colour</b>              | dark grey  |          |

### Indications

- Dorsolumbar pain following dorsolumbar vertebral compression fractures
- Dorsolumbar vertebral compression fractures of osteoporotic or metastatic origin
- Post-traumatic dorsal hyperkyphosis or due to arthritis
- Juvenile osteochondrosis
- Other indications as per medical prescription



ANATOMICAL  
FIT



LOW  
PROFILE



EASY  
ADJUSTING



CUSTOMIZABLE  
FIT



EASY  
PULLING



RADIO  
LUCENT



PATENT  
EP2502608



TELESCOPIC



SINGLE PATIENT  
MULTIPLE USE



PATENT  
PENDING



**Headquarter: ORTHOSERVICE AG**  
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland  
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89  
info@orthoservice.com · www.orthoservice.com



**ORTHOSERVICE DEUTSCHLAND GmbH**  
Flugstraße 8 · D-76532 Baden-Baden · Deutschland  
Tel. +49 (0) 7221 40 533 0 · Fax +49 (0) 7221 40 533 20  
info@orthoservice.de · www.orthoservice.de

**ORTHOSERVICE POLSKA Sp. z o.o.**  
ul. Warszawska 416a · 42-209 Częstochowa · Polska  
Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785  
info@orthoservice.pl · www.orthoservice.pl

**OSUK Ltd**  
Units 47&50 North East BIC  
Wearfield  
Sunderland Enterprise Park East  
Sunderland SR5 2TA  
Tel +44 (0) 191 51 66 220  
info@orthoservice.uk · www.orthoservice.uk

**RO+TEN S.r.l.**  
Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia  
Sede operativa e amministrativa:  
Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia  
Tel. 0039 0362 18 51 200  
info@roplusten.com · www.roplusten.com



orthoservice.com



**ORTHOSERVICE  
RO+TEN**

Take care feel better